


Lancing Device

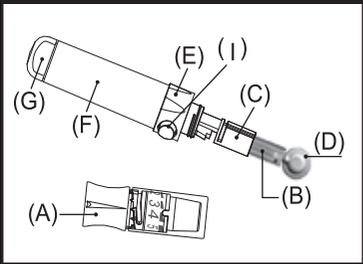

Instructions for Using Lancing Device

- | | | | |
|-----------|------------------------------|-----------|------------------------|
| A. | Adjustable end cap | E. | Device body |
| B. | Lancet | F. | Sliding barrel |
| C. | Lancet holder | G. | Color indicator |
| D. | Lancet protective cap | I. | Trigger button |

Precautions:




- Never share a lancet or a lancing device with anyone.
- Always use a new, sterile lancet-lancets are for single use only.
- Keep the lancing device and lancet away from latex, oil stains and dust.
- Do not immerse the lancing device in any liquid.

This lancing device is used to obtain capillary blood samples for blood glucose testing. The lancing device has 5-setting depth adjuster. Clean your hands thoroughly before using the lancing device. Use for self-testing.

Instructions for use:

1. Select a site on an outer edge of the finger or the forearm for the best blood drop. Wash the finger or the forearm with soap and warm water. Make sure the testing site is dried thoroughly.
2. Twist off the clear adjustable end cap (A). Insert the lancet (B) firmly into the lancet holder (C).
3. Remove the lancet protective cap (D) by holding the device body (E) and the lancet (B) firmly with one hand and twist off the lancet protective cap (D) with the other hand. (Note: Do not touch the lancet tip and save the lancet protective cap (D) for safe disposal of used lancet later.)
4. Replace the clear adjustable end cap (A) by putting it back onto the device body (E).
5. Set lancing depth by turning the arrow pointer to one of the five depth indicators on the end cap (A). For soft skin, please set the arrow pointer between 1-2 and for average skin, please select 3-4. For thicker skin, please select 4-5.
6. Select a soft, fleshy area on your forearm, away from bone and free of visible veins and hair. Hand your arm down and gently rub the selected area on forearm or from wrist to the base of finger. This will help improve blood flow.
7. Hold the end cap (A) in one hand and sliding barrel (F) in other hand. Gently pull back the sliding barrel (F) until a soft click is heard.
*Value-Added Feature: color indicator (G) in the transparent cap helps users to verify if the device is set ready for use.
8. Place the lancing device firmly against the selected finger site or the forearm, and press the round trigger button (I) to obtain blood sample. (Note: Gently massage the finger or the forearm before you obtain blood sample.)
9. After finishing the blood testing, remove the used lancet (B). Recap lancet using its saved lancet protective cap (D) and dispose the lancet in an appropriate container.

***To clean the lancing device and clear cap, wipe them with a soft cloth dampened with water and mild detergent. Do not immerse the lancing device in any liquid.**

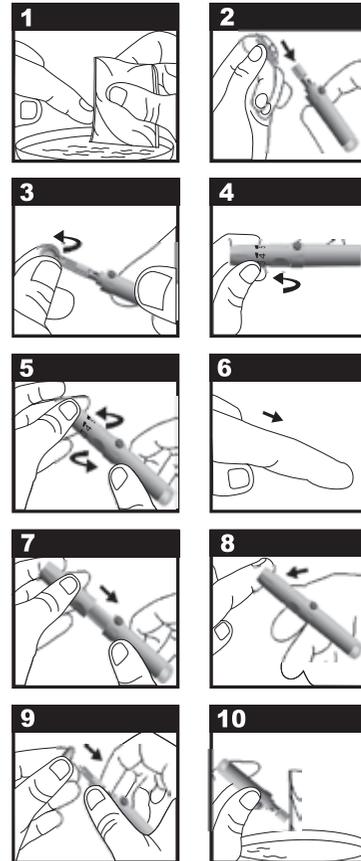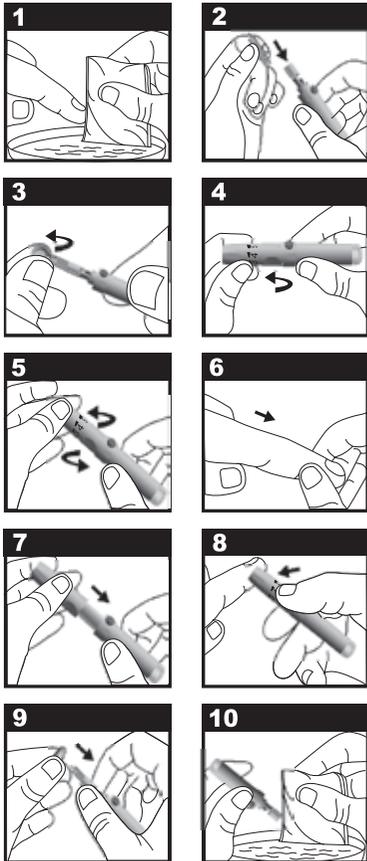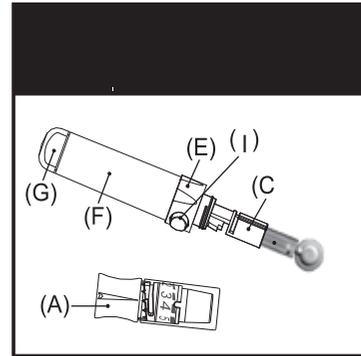

Distributed by:



Home Aide Diagnostics, Inc.

1072 S. Powerline Rd., Deerfield Beach, FL 33442, USA
TEL: 1-800-915-0116




Biotech Medical Corporation
New York, NY, USA



MEDNET GmbH
Borkstrasse 10, 48163 Münster
Germany

Rev.date:11/2012 8128HA-M02A2

Lancing Device

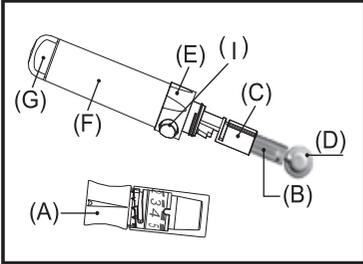

Instructions for Using Lancing Device

- | | | | |
|-----------|------------------------------|-----------|------------------------|
| A. | Adjustable end cap | E. | Device body |
| B. | Lancet | F. | Sliding barrel |
| C. | Lancet holder | G. | Color indicator |
| D. | Lancet protective cap | I. | Trigger button |

Precautions:




- Never share a lancet or a lancing device with anyone.
- Always use a new, sterile lancet-lancets are for single use only.
- Keep the lancing device and lancet away from latex, oil stains and dust.
- Do not immerse the lancing device in any liquid.

This lancing device is used to obtain capillary blood samples for blood glucose testing. The lancing device has 5-setting depth adjuster. Clean your hands thoroughly before using the lancing device. Use for self-testing.

Instructions for use:

1. Select a site on an outer edge of the finger or the forearm for the best blood drop. Wash the finger or the forearm with soap and warm water. Make sure the testing site is dried thoroughly.
2. Twist off the clear adjustable end cap (A). Insert the lancet (B) firmly into the lancet holder (C).
3. Remove the lancet protective cap (D) by holding the device body (E) and the lancet (B) firmly with one hand and twist off the lancet protective cap (D) with the other hand. (Note: Do not touch the lancet tip and save the lancet protective cap (D) for safe disposal of used lancet later.)
4. Replace the clear adjustable end cap (A) by putting it back onto the device body (E).
5. Set lancing depth by turning the arrow pointer to one of the five depth indicators on the end cap (A). For soft skin, please set the arrow pointer between 1-2 and for average skin, please select 3-4. For thicker skin, please select 4-5.
6. Select a soft, fleshy area on your forearm, away from bone and free of visible veins and hair. Hand your arm down and gently rub the selected area on forearm or from wrist to the base of finger. This will help improve blood flow.
7. Hold the end cap (A) in one hand and sliding barrel (F) in other hand. Gently pull back the sliding barrel (F) until a soft click is heard. *Value-Added Feature: color indicator (G) in the transparent cap helps users to verify if the device is set ready for use.
8. Place the lancing device firmly against the selected finger site or the forearm, and press the round trigger button (I) to obtain blood sample. (Note: Gently massage the finger or the forearm before you obtain blood sample.)
9. After finishing the blood testing, remove the used lancet (B). Recap lancet using its saved lancet protective cap (D) and dispose the lancet in an appropriate container.

***To clean the lancing device and clear cap, wipe them with a soft cloth dampened with water and mild detergent. Do not immerse the lancing device in any liquid.**

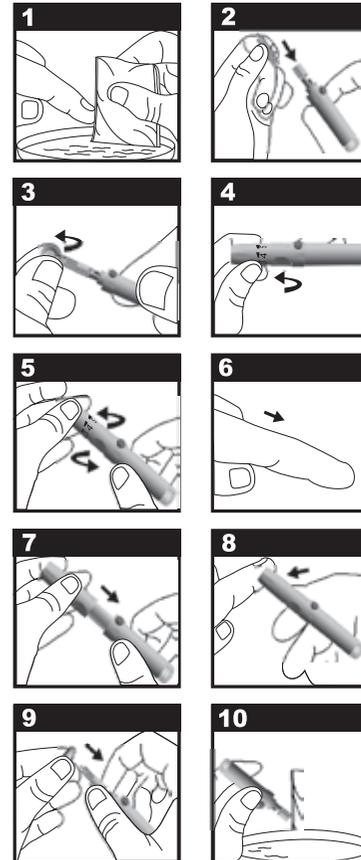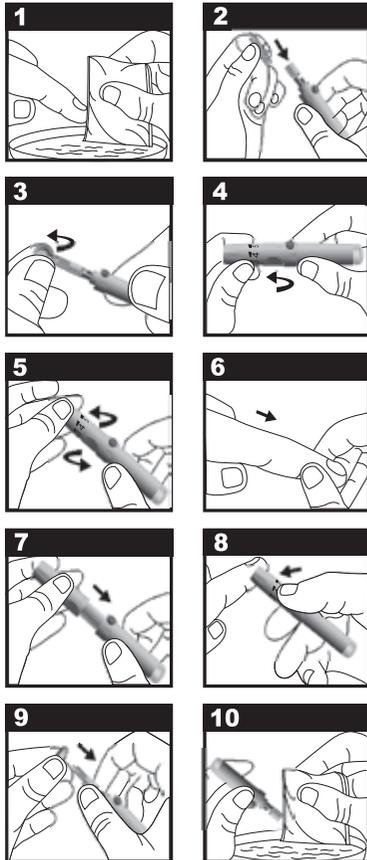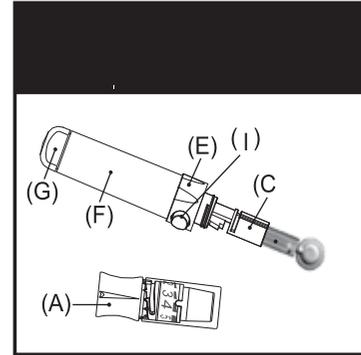

Distributed by:



Home Aide Diagnostics, Inc.

1072 S. Powerline Rd., Deerfield Beach, FL 33442, USA
TEL: 1-800-915-0116



Biotest Medical Corporation
New York, NY, USA
New Zealand, New Zealand
New Zealand, New Zealand



MEDNET GmbH
Borkstrasse 10, 48163 Münster
Germany

Rev.date:11/2012 8128HA-M02A2